

Electronic Balances

TW/TX/TXB series

The beginning of the new standard.
Extremely capable, but easy to operate.



Internal Calibration

Press just two keys to calibrate the balance whenever calibration is necessary (TW only). Calibration is very fast, taking only 15 sec.




Easy Setting

Quickly adjust the desired ratio of stability and response for every application, even during measurement, with one-touch operation.



Menu Operation Key

Menu navigation keys are separated from weighing operation keys and arranged in a familiar 5-way navigation circle. Up, Down, Right, Left and Enter are the simple operational steps.

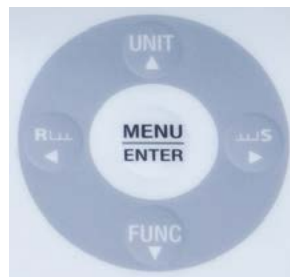

Dry Battery Operation

Power the TXB series balances with an AC adapter or batteries.



Piece Counting

Unit weights of up to five different samples can be easily entered, stored and recalled for use.



Measuring Carats

In addition to grams (g), weigh in %, number of pieces, ct, kg, mg, lb, oz, TTI, etc. or a custom conversion unit, more than 20 units in all. Change quickly from display of % or counting to gram weight (g) display.



Checkweighing

Compare samples to target values or pass/fail criteria and clearly indicate the results.

TW Series



TX Series



TXB Series



TW223L TW323L TW423L
TX223L TX323L TX423L

TW Series

Model name	TW223L	TW323L	TW423L	TWC323L	TWC623L
Capacity	220 g	320 g	420 g	320 ct (64 g)	620 ct (124 g)
Minimum Display	0.001 g	0.001 g	0.001 g	0.001 ct (0.0002 g)	
Pan Size (mm)	Approx. ø110			Approx. ø80	
Dimensions	Approx. W206 × D291 × H241 mm				
Weight	Approx. 4.2 kg			Approx. 4.1 kg	



TX2202L TX3202L
TX4202L

TX Series

Model name	TX223L	TX323L	TX423L	TX2202L	TX3202L	TX4202L	TXC323L	TXC623L
Capacity	220 g	320 g	420 g	2200 g	3200 g	4200 g	320 ct (64 g)	620 ct (124 g)
Minimum Display	0.001 g	0.001 g	0.001 g	0.01 g	0.01 g	0.01 g	0.001 ct (0.0002 g)	
Pan Size (mm)	Approx. ø110			Approx. W167 × D181			Approx. ø80	
Dimensions	Approx. W206 × D291 × H241 mm			Approx. W200 × D291 × H80 mm			Approx. W206 × D291 × H241 mm	
Weight	Approx. 3.8 kg			Approx. 2.8 kg			Approx. 3.8 kg	



TWC323L TXC323L
TWC623L TXC623L

TXB Series

Model name	TXB422L	TXB622L	TXB4201L	TXB6201L
Capacity	420 g	620 g	4200 g	6200 g
Minimum Display	0.01 g	0.01 g	0.1 g	0.1 g
Pan Size (mm)	ø110		ø160	
Dimensions	Approx. W199 × D260 × H77 mm			
Weight	Approx. 1.5 kg			



TXB2201L TXB6201L
TXB4201L TXB6200L



TXB222L TXB622L
TXB422L TXB621L

Optional Accessories

Description
Electronic Printer EP-100
Electronic Printer EP-110
RS-232C Cable
In-use Protective Cover
In-use Protective Cover for Display
USB Conversion Kit



Electronic Printer
EP-100



Electronic Printer
EP-110